



SIDDHI SHAKSHI

PHARMA EQUIPMENT



We are Manufacturer of,
**Pharmaceuticals, Chemicals, Cosmetics,
& Food Processing Machineries**

E-mail : info@siddhishakshipharma.com

Web: www.siddhishakshipharma.com

Contact : Mr. Anil Sharma

Cell: +91 7875328690, 9028808609

MSME
Registered



Fluid Bed Dryer (FBD)

A fluidized bed dryer (also called fluid bed dryer) is a kind of equipment used extensively in the pharmaceutical industries to reduce pharmaceutical powder and granules moisture content. Drying is essential for wet granules for compression into the tablet and modifying viscous and socky materials' characteristics. Drying is commonly the last step in the unit process before compression into tablet form and packaging.

Salient Features :

- Very Versatile Process Equipment for Uniform and Efficient Drying of Products.
- Fully Automatic Operation Based on PLC Control / Without PLC Control
- Minimum Handling and more Hygienic
- Lifting & Tilting Device of Product Container for Easy Handling
- Capacity Available from 5 Kgs to 500 Kgs.
- FLP Models available as Option
- MOC S.S. 304/ 316/ 316L
- Option – Steam Heated or Electric Heated Blower
- CGMP Model with Documentation (IQ, DQ, PQ, OQ) – As Option.



Octagonal Blender

Salient Features :

- Suitable for Quick and Efficient Mixing for a Complete Range of Powders, Particulate and Granular Materials.
- Homogeneous Mix. Delicate Tumbling Action with a Specially Designed Baffle.
- Easy Charging of Materials through S.S. Canisters.
- Capacity – 20 Ltrs. to 4000 Ltrs.
- MOC – S.S. 304/316/316 L
- Shape of Container – Polygon Shaped Frustum.
- S.S. Octagon type body with Rectangular Central Portion Top and Circular Bottom with a Butterfly Valve.



SIDDHI SHAKSHI PHARMA EQUIPMENT

For More Please Contact Us:

E-mail : info@siddhishakshipharma.com

Web: www.siddhishakshipharma.com

Contact : Mr. Anil Sharma

Cell: +91 7875328690, 9028808609

MSME
Registered



Double Cone Blender

Salient Features

- The conical shape at both ends enables uniform mixing and easy discharge.
- The cone is statically balanced which protects the gear box and motor from any excessive load.
- Powder is loaded into the cone through a wide opening and discharged through a butterfly or a Slide valve.
- Depending upon the characteristic of the product, paddle type baffles can be provided on the shaft for better.
- Contact parts are made of SS 304 or SS 316.
- 'Slant' design (off centre) CLIN CONE BLENDER.
- Dust free bin charging system ensures minimum material handling.



Starch Paste Kettle

Salient Features

- It is Specially Designed for Preparation of Binder Material.
- Hemispherical Kettle with Jacket Mounted on Trunnion Supports with the Help of Shaft.
- The Kettle has a Small Cylindrical top Shell With a Spout for Pouring out Binder Mass.
- Jacket is Provided with Steam Accessories, Pressure Gauge, Safety Valve and Vent Cock.
- The Kettle is also Provided with a Lever with Hand Wheel to Accomplish Tilting Manually.
- The Kettle is Provided with 16 Swg Thick
- S.S.304 Loses Lid.



Coating Pan



Salient Features :

- Facility of angular adjustment of coating pan facilitate the perfect angle for coating deferent product.
- Available Sizes 12" 24" 30" 36" 48" 60" 72"

SIDDHI SHAKSHI PHARMA EQUIPMENT

For More Please Contact Us:

E-mail : info@siddhishakshipharma.com

Web: www.siddhishakshipharma.com

Contact : Mr. Anil Sharma

Cell: +91 7875328690, 9028808609

MSME
Registered



Vibro Sifter

Salient Features

- MOC : SS316L / SS316 / SS304
- Flame proof motor with FLP push button switch.
- Magnetic arrangement for trays of iron particles.
- Extra deck available
- GMP Model



Multimill

Salient Features :

- Multi speed drive for different item pounding.
- Simple destroying and cleaning.
- High level of surface completion.
- Simple convenience.
- MOC SS 304/316/316L.
- cGMP with documentation



Rotocone Vacuum Dryer



Rotocone Vacuum Dryer is suitable for drying materials that cannot resist high temperature, easily oxidized materials, volatile materials that should be retrieved, and potent irritants and poisonous in nature. The Rotocone Vacuum Dryer with improved technology integrates during operation under vacuum.

The Rotocone Vacuum Dryer facilitates enhanced drying efficiency, low-temperature operation, and process economy by total solvent recovery. It helps cGMP base working by achieving optimum dust control while offering efficient charging and discharging of materials. The drying unit equipped with lump breakers initially breaks large lumps and subsequently powders them (in a single cone). The dryer's rotary action, together with the breakers' mechanical action, cuts down drying time and gives a lump-free product.

SIDDHI SHAKSHI PHARMA EQUIPMENT

For More Please Contact Us:

E-mail : info@siddhishakshipharma.com

Web: www.siddhishakshipharma.com

Contact : Mr. Anil Sharma

Cell: +91 7875328690, 9028808609

MSME
Registered



Ribbon Blender

We have become the leading association involved in offering an extensive array of Talcum Powder Manufacturing Plant / Talcum Powder Mixing Machine to our honored clients. It is manufactured using an excellent quality mechanism and pioneered technology as per the set industry standards. This machine is carefully observed on numerous parameters by our quality experts to eradicate any flaw.

Keeping the excellence in this field, We manufacture and supply the Talcum Powder Mixing Machine's best quality. Due to its excellent mixing capacity and robust construction, the offered powder mixing machine is highly demanded in Indian and International markets. Our customers can avail of this mixing machine in various technical specifications.


Storage Vessel

We have been supplying stainless steel tanks and process equipment to the brewing, beverage, food, waste, pharmaceutical and chemical industries. We also supply mild steel tanks.

We handle the increasing demand for high quality stainless steel storage and process vessels from a wide range of industries - each with its own requirements.



Syrup Manufacturing Plant

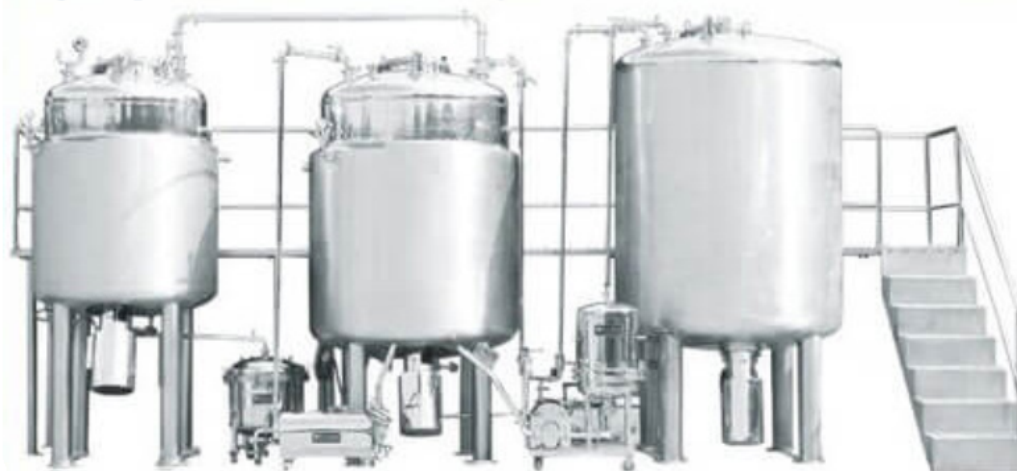

Salient Features :

- CAPACITY : 50ltr - 7000 ltr
- Sugar and Glucose Transfer System
- Sugar Syrup Preparation with proficient base passage unsettling.
- Clear and Suspended Liquid Manufacturing System with On-line or Built-in Homogenizer
- Proficient Filtration System for the unmistakable/suspension syrups
- CIP and SIP System for cleaning and purification of the framework

SIDDHI SHAKSHI PHARMA EQUIPMENT

For More Please Contact Us:

E-mail : info@siddhishakshipharma.com

Web: www.siddhishakshipharma.com

Contact : Mr. Anil Sharma

Cell: +91 7875328690, 9028808609

MSME
Registered



Zero Holdup Sparkler Filter Press

Salient Features :

- The complete line is manufactured under cGMP standard.
- All contact parts are made of ss316.
- The plant is designed to be operated by one operator and one helper saving on precious manpower cost.
- Available sizes in 8-33 plates



Pressure Nutsche Filter

Salient Features :

- 316 Stainless Steel fabrication or other alloys including Hastelloy C-276/C-22 optionally available.
- For aqueous regia production nutsche filter available with halar or tefzal coating.



Pressure Vessel

Pressure vessels are leak-proof containers that store liquid or gas. Pressure vessels of various sizes and shapes have been produced for different purposes. Generally, preferred geometries are spherical, conical, and cylindrical. A typical model is the combination of a long cylinder with two heads. Pressure vessels work at internal pressures higher or lower than air pressure.



SIDDHI SHAKSHI PHARMA EQUIPMENT

For More Please Contact Us:

E-mail : info@siddhishakshipharma.com

Web: www.siddhishakshipharma.com

Contact : Mr. Anil Sharma

Cell: +91 7875328690, 9028808609

MSME
Registered